

V8v Electric bed


Registration code: SKD-C



SAIKANG/Medical



V8v FEATURES AND HIGHLIGHTS

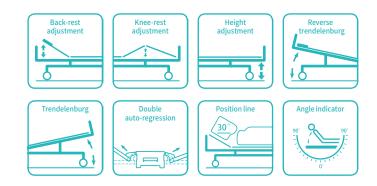

Technical parameters

• External size (LxW):	2085x1025 mm
Mattress platform:	1925x900 mm
Safe working load:	250 kg
Height adjustment:	455-745 mm
Back-rest adjustment:	0-75° (土5°)
Knee-rest adjustment:	0-32°(土5°)
Reverse trendelenburg:	0-12°
Trendelenburg:	0-12°
Angle indicator:	0-90°
Auto-regression:	67+57 mm
Position line:	30°
 Voltage and frequency: 	100-240 V · 50/60 Hz
 Plug as required (hospital grade) 	

Technical configuration

🗹 Linear motor	4 pcs
☑ Hand controller	1 pc
☑ 5'' double side castors, central locking	4 pcs
☑ PP side rail extra height	1 set
\blacksquare PP bed ends with safe lock with Φ 65 bumper wheels	1 set
PP bed platform with integrated mattress retainer	1 set
\blacksquare Patent angle indicator for back-rest and trendelenburg position	4 pcs
☑ IV pole 4 hooks	1 set
☑ IV pole prevision	4 pcs
🗹 Drainage hook	2 pcs
🗹 Central locking pedal	2 pcs

Optional

Emergency stop button	
Trapeze bracket/pole	

Optional color (MOQ100pcs- color changing)



POSTURE FUNCTION



Safe working load

 250_{kg}



Height adjustment





Back-rest adjustment

SAIKANG

75°



Reverse trendelenburg





Knee-rest adjustment

32°



Trendelenburg

QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating. efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



Design

The unique germ-free design meets hospitals' strict infection control requirment. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.



Plastic | Painting Environmental plastic can be sterilized in 100°C, tensile upto 30MPa. 11 Process epoxy painting, ASTM testing anti-baterial, paint thickness 0.12mm, brightness 60°, paint can seric t 50kr impact



Auto-Regression Auto-Regression 10+9 cm When repositioning the backrest and thigh rest, the V bed makes the back plate retreat 10cm and the thigh plate retreats 9 cm so to minimises the pressure and shear forces in the patients' pelvic region. This results in a significant reduction in the development of pressure ulcers.



Steels | welding Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



Electric Motor American motor, with UL/ROHS/EN standard. Noise under 50dB, with

lifetime 20,000 times movement.



Bed Frame Solid bed base made of 5.0mm stamped rotating plate.

Accessories



Platform

Integrated PP bed platform, adopting advanced manufacturing technology. 4-5 section platform, enables independent installation and easy maintenance.

Enlarge the seat area when back part is raised, providing more comforts for the patients.

Reduce the depression of the sponge at the seat panel.



Castors

 $5^{\prime\prime}$ Double side castors, central locking. TPR tire no worn out after running 10 km.

Anti-winding hard shell cover, united forming without bolts.

 Pass dynamic test: pass obstacles 40000 times, under the standard 200kg load of the whole bed.

Head & Foot Board

LG Chem fresh PP material, Plug in type. Equipped with the bumpers wheel, diameter 65mm, preventing collision.







Tuck away side rail with buffer integrated. Integrated handle facilitates the patients to get up/off the bed.





Wall bumper wheel

Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall.



Hand Controller Hand controller: flexible spiral cable for easy mobility and operation in optimal position. The nurse can control all standard functions.





Platform integrated mattress retainer at the both sides of the back and leg section, with a boarder of 33mm height.













Design





Manufacturing



Maintenance



Service



South America Branch

Saikang Medical Technology, Inc W3400 Corrientes, Argentina TEL: +55 (11) 2758-0438 Rua Universal, Hollywood, São Bernardo do Campo, Brazil TEL: +55 (11) 2758-0438

East Africa Branch Saikang Medical Technology Ltd Mezzanine Eloor Tai Towe

Mezzanine Floor, Taj Tower, upper hill, nairobi, Kenya TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd Sala Thammasop Subdistrict Thawi Watthana District Bangkok Thailand TEL: +66 642741046

US Branch

Healthward Technology, INC 17800 CASTLETON ST STE 180. CITY OF INDUSTRY CA 91748 TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG' 454004, Russian Federation, Chelyabinsk, Universitetskaya Naberezhnaya, Russia TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC Musaffah P.O Box 53200, Abu Dhabi, UAE Tel: +971 585710652



Jiangsu Saikang Medical Equipment Co., Ltd. STOCK CODE: 870098

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China +86 18021231901 · S 24/7 export@saikangmedical.com

www.saikangmedical.com





o General catalogue